



the difference is clear

Peak Flow Meter asmaPLAN & Peak Flow Meter asmaPLAN+

The simple, accurate way to monitor and manage asthma

PEF asmaPLAN

- Accurate
- Compact, lightweight and portable
- Robust and hard wearing
- Features comfortable integral mouthpiece
- Easy to read scale markers
- For adults and children
- Multi-subject use with SafeTway[®] mouthpiece
- Conforms to all peak flow standards including EN13826:2003

PFM asmaPLAN+


Once a patient's target peak expiratory flow (PEF) has been established by a healthcare professional the color coded dual sliders on the front of the asmaplan+ can be set. The subject personal action plan can be written on the coordinating action plan label on the back-side make self-monitoring and self-management simple for asthma patients.

Vitalograph is the world's leading provider of outstanding respiratory diagnostic, screening and monitoring devices. With a heritage of excellence spanning half a century, Vitalograph continue to make pioneering contributions to effective respiratory care and enhanced quality of life.





the difference is clear

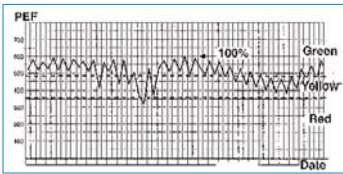

Technical Specification:

Type of Device
Variable orifice simple mechanical
peak flow meter

Peak Flow Meter Standard:
EN13826:2003

Safety Standards:
EN ISO 60601-1:2005 -1 -2

Medical Safety Standard
Medical Devices Directive
93/42/EEC 1993

Designed & Manufactured Under:
ISO 13485:2003, FDA 21CFR820,
CMDR & JPAL

Range: (L/min):
50-800L/min

Accuracy:
Better than 10%

Repeatability:
Better than 5% or 10L/min whichever
is greater

Inter-Instrument Variation:
Better than 5% or 10L/min whichever
is greater

Scale Intervals:
10L/min below 700, 20L/min above

Daily Record Chart
Included

Measuring Principle:
Piston/spring and moveable pointer
Frequency Response:
A/B difference less than 15L/min/15%
(Annex B of EN13826)

User Instructions:
Multi-language instruction sheet

Clinic Use:
Certify accuracy at least annually
using a certified calibrated
spirometer.

Home Use:
Replace every three years

Net Weight
2.6 oz

Highest Resistance to Flow:
0.00384KPa/L/min@ 720L/min

Storage Temperature
32-122°F

Relative Humidity: 10%-95%

Dimensions:
6" x 2 1/2" x 1 3/8" (including integral
mouthpiece)

Ordering Info:

Catalogue No.

43600	Peak Flow Meter asmaPLAN
43700	Child Peak Flow Meter
43800	Peak Flow Meter asmaPLAN+
20242	SafeTway® mouthpieces (200)
20980	Mini SafeTway® mouthpieces (50)

Vitalograph Inc.
13310 West 99th Street
Lenexa
Kansas 66215 USA
Tel.: (913) 888-4221
Fax: (913) 888-4259
E-mail: vitcs@vitalograph.com

Vitalograph GmbH.
Jacobsenweg 12
22525 Hamburg
Germany
Tel.: (040) 547391-0
Fax: (040) 547391-40
E-mail: info@vitalograph.de

Vitalograph Ltd.
Maids Moreton
Buckingham
MK18 1SW England
Tel.: +44 (0) 1280 827110
Fax: +44 (0) 1280 823302
E-mail: sales@vitalograph.co.uk

Vitalograph (Ireland) Ltd.
Gort Road Business Park
Ennis
Co. Clare Ireland
Tel.: (065) 6864100
Fax: (065) 6829289
E-mail: sales@vitalograph.ie



www.vitalograph.com

Vitalograph®, Spirotrac®, SafeTway®, BVF™, asmaPLAN™ are trademarks of Vitalograph



PRINT REF: 14429_1